

Manual

Version 20.03.2023

Transmitter Model C1

Turbidity Probes Models AP2, TURBImat-HC, TURBImat-LC,
Absorptionphotometer

Chemtronic Waltemode GmbH

40789 Monheim am Rhein
Niederstr. 14
Germany

Phone: +49-(0)2173-57007 Fax: +49-(0)2173-56007
E-Mail: info@chemtronic-gmbh.de
Internet: www.chemtronic-gmbh.de



Content

Copyright.....	3
Non-warranty.....	3
Warranty conditions	3
Storage.....	3
Safety instructions	4
For ALL chemical liquids	4
Electrical installation.....	4
Hazardous area.....	4
Maintenance	4
Package Content.....	4
Installation Guidelines.....	5
Ideal Installation for Probes Models AP2 and TURBImat-LC	6
Ideal Installation for Probe Model TURBImat-HC	6
Principle of Operation	7
What does Turbidity mean?.....	7
What Causes Turbidity?	7
Measurement of Turbidity.....	7
Principle of Absorbance Measurement	7
Absorption Photometry	8
Turbidity Probe Model TURBImat-HC	9
Software Installation (C1-Config Software)	10
Connect C1 Transmitter to a PC- or Laptop	10
C1-Config Software	11
Select Language	11
Select COM Port	12
Set measuring Range.....	13
Enter measuring Unit	13
Set Damping	13
Adjust Zero Offset	13
Adjust Gain	13
Set Relay Setpoint.....	14
Set Relay Hysteresis	14
Apply Settings to Transmitter	14
Save Settings to File.....	14
Read Settings from File.....	14
Calibration.....	15
Measuring Screen	17
Maintenance Work	18
List of Spare Parts / Wearing Parts.....	18
Trouble Shooting	19
Terminal PC-Board Model C1	20
Wiring Diagrams Model AP2, TURBImat / C1.....	21
Dimensional Drawing Model C1	22
Dimensional Drawings Model AP2.....	23
Dimensional Drawings Model AP2 with Flange.....	24
Dimensional Drawings Model TURBImat-HC.....	25
Dimensional Drawings Model TURBImat-LC	27
Data Sheet Transmitter Model C1	28
Data Sheet Probe Model AP2.....	29
Data Sheet Probe Model TURBImat-HC.....	30
Data Sheet Probe Model TURBImat-LC.....	31
EU Declaration of Conformity (CE)	32

Copyright

© 2023 Chemtronic Waltemode GmbH

This manual including all of its parts are protected by copyright of Chemtronic Waltemode GmbH.

Any further use beyond the copyright laws is not allowed without the written approval of

Chemtronic Waltemode GmbH.

There are no further reaching warranty claims to take due to the content of existing document.

Non-warranty

Chemtronic Waltemode GmbH takes no responsibility or warranty for technical or editorial mistakes or omissions in this document. The information in this publication will be provided without warranty for their correctness. Chemtronic Waltemode GmbH reserves all rights to make changes to this document without announcement.

Chemtronic Waltemode GmbH is by no means liable for any consequential damages or miscellaneous damages, which occur during the use or storage of the system.

Some countries do not allow the limitation of liability, so that this exclusion is possibly not valid for you.

Warranty conditions

The warranty for this product takes place according to the legal provisions. As place of warranty we explicitly name Chemtronic Waltemode GmbH, Monheim am Rhein, Germany. Chemtronic will charge the additional expenses e.g. travel expenses, overnight cost, etc., in case of warranty service at the customer site.

Consumable parts as well as wearing parts e.g. lamps, gaskets, etc. are not covered by warranty.

Improper use, gross negligence, use of the system outside of the technical specifications, missing or improper maintenance will cause a loss of the warranty claim. Please pay attention to our general terms of contract.

Storage

Please inspect the instrument immediately after receiving for potential shipping damages.

If the instrument has already been unpacked for inspection or testing, or if the instrument has been removed from the process and it is not to be installed or reinstalled for more than 1 day, the following procedure should be observed:

1. If the instrument has been in service, the wetted portion should be thoroughly cleaned (typically with clean water) and then thoroughly dried.
2. In case the original packing material is not available, place the instrument in a sealed heavy plastic bag with a desiccant added to assure clean dry storage.
3. The instrument should then be stored in a protected area until time of installation.

Safety instructions

Take notice to the following general safety instructions during use and operation of the system. Ignoring these instructions or special warnings inside of this manual can damage the sensor, cause inaccurate measurements, and possibly result in unsafe installations. Chemtronic Waltemode GmbH will not take any responsibility for consequences arising from ignoring the safety instructions and warnings.

For ALL chemical liquids

Reagents, standards or the measured product may be toxic or otherwise harmful to human health. Observe appropriate precautions when handling and disposing the reagents. Review the appropriate published material safety data sheets regarding the safe handling and disposal of all these products.

Electrical installation

Qualified technical personnel must perform the electrical installation of the system.

Hazardous area

DO NOT INSTALL the system in hazardous.



Do not install the system in hazardous areas (Ex Zone 0, Ex Zone I, Ex Zone II) without the optional explosion proof approvals and certifications.



Maintenance

Always disconnect the instrument from power during maintenance, replacement of components, installation of additional components or any other operations at the open instrument.

Only qualified technical personnel must perform this work.

Operating the instrument with open enclosure

Only qualified technical personnel should operate the instrument when the enclosure is open (e.g. during repair or calibration procedure). Be careful that no moisture enters the enclosure.



Some components inside the instrument are energized with voltages, which can cause lethal shocks in case of contact. Be careful during installation, handling and operation of the instrument.



Improper installation / operation of the system

Warranty is void if the system is / was installed improperly, handled improperly, used outside of the technical specifications or damaged by gross negligence.

Package Content

- 1.) Transmitter model C1 with Turbidity probe according to users specification.
- 2.) USB Stick with TURBIsmart-F Config- / Calibration Software / Manuals and Certs
- 3.) USB cable.
- 4.) Printed Manual and Documentation



Please contact us immediately in case of an incomplete delivery or in case of an damage off one of the components.



Installation Guidelines



The measuring system is typically delivered with a factory calibration according to user's specification. The device is ready for use without further calibration / configuration.



Make sure that Instrument is disconnected from the mains, during wiring procedure.

Danger, take notice!



Make sure that your pipe is not pressurized or at low pressure before you insert or retract your probe to process. Before beginning any work the sensor must be cleaned and flushed thoroughly. Depending on customers' application product residues can be very dangerous (aggressive, poisonous).

Wear protective equipment if necessary (masks, safety glasses, protective gloves, etc.).

Transmitter Model C1

- Transmitter and Sensor are sensitive components, do not drop and handle with care.
- Please do not use the system in case of any visible mechanical damages!
- Do not operate with open enclosure. Exposed electronics will promote the risk of electrical shock and might cause damage to the unit.
- The enclosures of sensor and transmitter fit IP65 / Nema 4x protection class.
- Follow the wiring diagram in this manual.
- Supply voltage: 24VDC+/-10% (optional 100-240VAC, 50-60Hz), observe type plate!
- Proper instrument grounding: The instrument must be grounded to avoid electrical shock hazards to personnel and damages to the unit.
- The power consumption of the Instrument is less than 15W.
- Relay ratings for alarms and control of external devices are 1 A at 24VDC.
- Analog output: 0/4-20mA, max. 500Ohm (isolated)
- The ambient temperature for the Transmitter must be in a range of -10°C to 60°C (use dry purge air above 40°C approx. 10 l/h).


Probe Model AP2, TURBImat-HC or TURBImat LC

- The instrument is typically factory pre- calibrated according to customers specification please double check before installation!
- Avoid air and gas bubbles in your process, they cause measuring disturbances. Air and gas bubbles will cause noise and drift of the measuring signals. (The air bubbles are not expected at pressures upwards of 2bar in aqueous solutions).
- The process pressure of pipe insertion probes respectively the immersion depth (open channel, basin or tank installation) should never exceed the specification of the sensor.
- The process temperature should never exceed specification of the probe.




Exceeding the specified maximum pressure and/or the specified maximum temperature represents a high safety risk.

Ideal Installation for Probes Models AP2 and TURBImat-LC

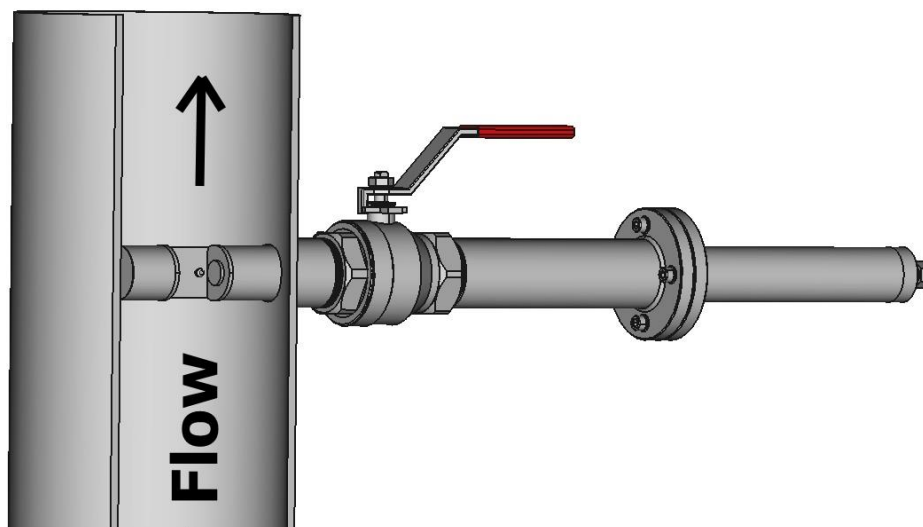

Ideal Installation for Probe Model TURBImat-HC

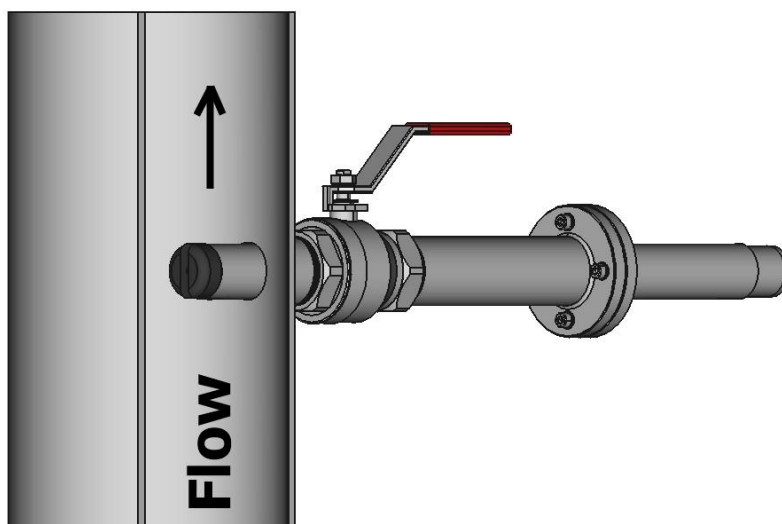

Principle of Operation

Absorption photometers will be used to measure the turbidity of liquid products.

The light absorption is measured at a specific wavelength, typically in the near infrared (NIR) at about 700nm to 850nm. The absorbance values are used to calculate the turbidity in an application specific measuring unit.

What does Turbidity mean?

Turbidity is an optical impression, which describes the characteristic of a transparent product, to scatter and absorb light. A focused light beam will be attenuated and scattered in hazy products, so that this product can become practically opaque in bigger layers.

What Causes Turbidity?

Turbidity is caused by particles in transparent products. A particle is defined as something with a different refractive index as the carrier product. Some examples of particles are minerals, yeast cells, metals, oil drops in water, milk in water, gas bubbles and aerosol's.

Measurement of Turbidity

Turbidity is not a clearly defined magnitude like e.g. temperature or pressure.


Turbidity is a subjective impression. For this reason turbidity measurement systems will be typically calibrated by using a comparison's standard such as Formazin or Diatomaceous Earth.

Principle of Absorbance Measurement

A light source and a detector are located parallel over against. The changing of light intensity, caused by the product between lamp and detector, will be detected and displayed as measuring result. The application is the measurement of concentrations (color or solids) in a liquid. The basic measurement unit of an absorption photometer is called AU (Absorption - Unit).

Definition of AU:

- 1 AU = 90 % loss of light intensity
- 2 AU = 99 % loss of light intensity
- 3 AU = 99,9 % loss of light intensity
- 4 AU = 99,99 % loss of light intensity
- 5 AU = 99,999 % loss of light intensity

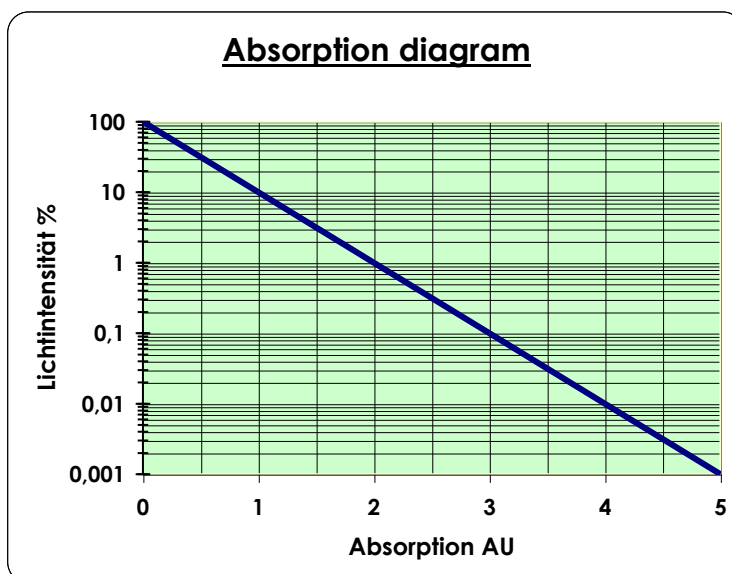

Typical Measuring units

ppm:	Parts per million	FNU ¹ :	Formazin Nephelometric Unit
FTU:	Formazin Turbidity Unit	mg/l:	Milligram per liter
TEF:	Trübungseinheiten Formazin (German for FTU)	gr/l:	Gram per liter
EBC:	European Brewery Convention	% TS:	Percent total solids
NTU ¹ :	Nephelometric Turbidity Unit		

The dependencies on the different measuring units

1 FTU = 1 TEF = 1 NTU¹ = 1 FNU¹ = 0,25 EBC

¹ Nephelometry describes the method of side scatter turbidity measurements, these units are used at 90° side scatter turbidity meters only.

Based on comparisons measurements, by using a 12° forward measurement system we have found the following dependencies.

1 FTU = 1 TEF = 0,25 EBC = 2,05ppm = 2,05 mg/l = 0,00205 g/l = 0,00000205 % TS

* At a specific particle weight of 1 kg/dm, 1mg/l particles in 1 kg of water will correspond to 1ppm.

Typical ranges

The original design of scatter light turbidity meter was used for the detection of low turbidities. The resolution of these instruments is suited easily in ranges lower as 0.1ppm (approx. 0.05 TEF / FTU / FNU / NTU or approx. 0.01EBC) and better. The maximum range is in ideal case lower than 200ppm, but there are as well systems available with a range of more than 8000ppm.

Please take notice:

There are turbidity meters with several different methods of measurement available!

Make sure that your turbidity meters operate according to the same principle of measurement, if you compare measurement results.

For example:

you can compare 90° vs. 90° scattered light turbidity meters, or 12° vs. 12° scattered light turbidity meters.

Never compare turbidity meters with different methods of measurement!

Absorption Photometry

The sensitivity of an absorption turbidity measurement system, is mainly dependent by two parameters.

1. The intensity of the light source, which is typically constant.
2. The optical path length (OPL¹): The optical path length is variable and can be changed to match the requirements of specific applications.


To measure low concentrations, a large OPL is required. Otherwise there is not enough change in light intensity, to get good measurement results.

To measure high concentrations, a small OPL is required. Otherwise there is not enough light energy, to pass the product.

¹ OPL [= optical path length] specifies the product slot between light source and detector


Turbidity Probes Models AP2 / TURBImat-LC

The probe model AP2 uses the single beam absorbance measurement to detect turbidity in liquids: The measurement wavelength is typically 880nm (other wavelength on request).

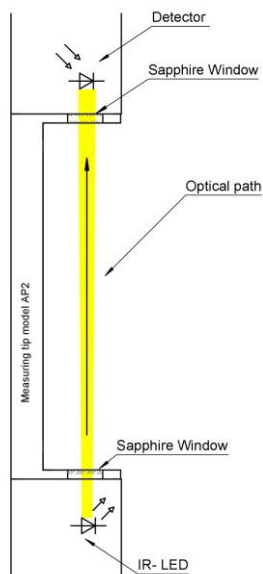

The beam of the IR- LED passes the Sapphire windows embedded in the sensing head. The detector captures the IR- absorption of the liquid in the optical path.

Turbidity Probe Model TURBImat-HC

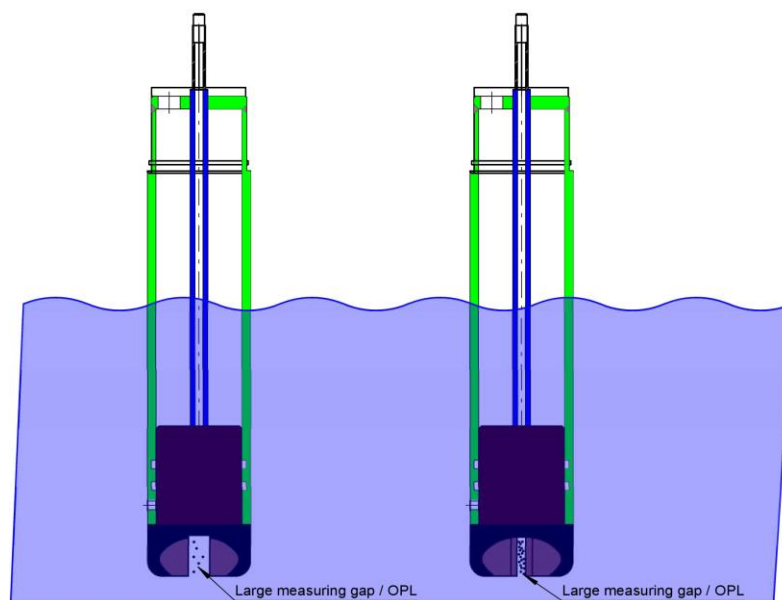

Large OPL = Low concentrations / high sensitivity

Small OPL = High concentrations / low sensitivity

Software Installation (C1-Config Software)

System Requirements

Supported Operating Systems	Microsoft Windows 7, 8, 10
Requirements	Microsoft.NET Framework 4.72 (x86 & x64)
RAM	2GB
Hard Disk Space	< 1MB

Installation Procedure C1-Config Software

- I. Close any open programs that are running on your computer before starting the installation.
- II. Double click on setup.exe.
- III. Software installs automatically

Connect C1 Transmitter to a PC- or Laptop

Connect the device according to the connection diagram and switch on the device. Then connect the PC/laptop to the C1 measuring amplifier using the USB cable (not vice versa).

To perform a restart, always disconnect the USB plug as well.



Take Notice:

The USB- connector of the C1 is located under the cover cap at the outside of transmitter enclosure.

C1-Config Software

Start Software per click/double click on the C1-Config icon.

Select Language




Chemtronic Waltemode GmbH
Niederstr. 14
40789 Monheim am Rhein
Germany


C1- Configuration Software



Select your language and proceed!

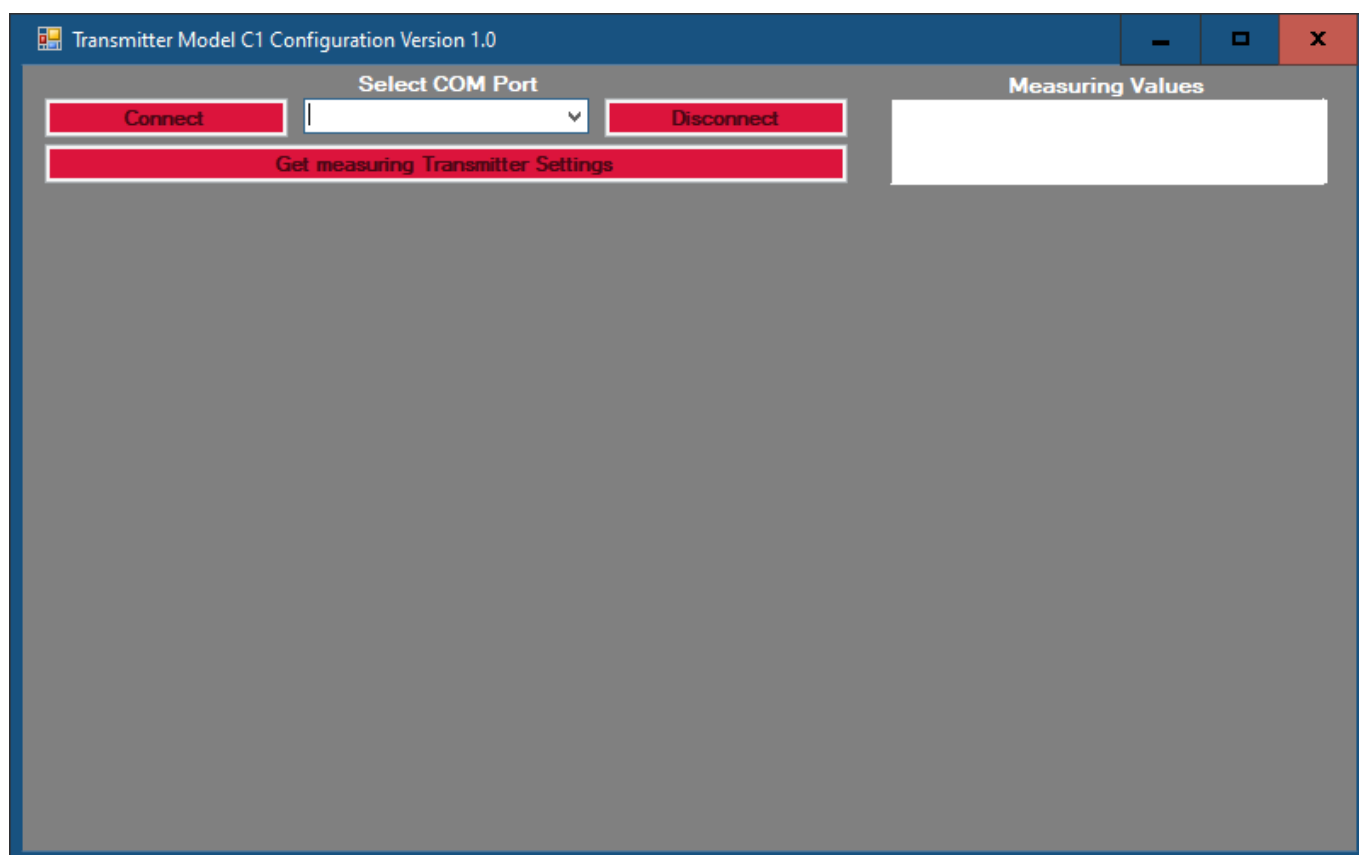

Select COM Port

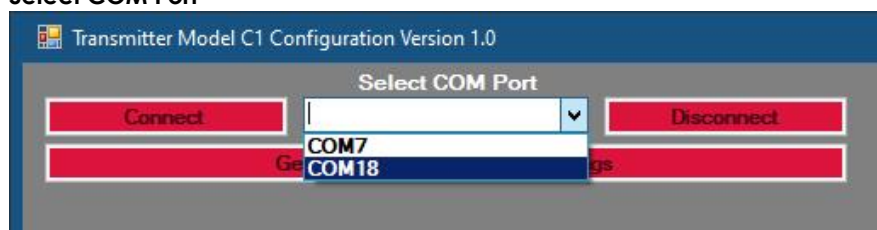

Select the COM port of your turbidity meter, click “connect” button and then the “Get Settings” button.



Finding the correct COM Port:

Shut down C1-Config software. Disconnect USB cable.

Restart software, go to “Select COM Port”, look to the list of active Ports.

Shut down C1-Config software again. Reconnect USB cable.

Restart software; go to “Select COM Port”, Com Port list shows now the additional Port of the instrument. Select these COM port and click “Connect Button”.

A click to the “Get measuring Transmitter Settings” button will transfer the current transmitter settings and measuring values to the screen.

Set measuring Range

Adjust or enter your measuring range.

Enter measuring Unit

Allows entering an application specific measuring unit like ppm, mg/l, EBC, FTU, ...
The measuring unit is limited to four characters (ASC II).

Set Damping

The Damping function averages the last measuring results of the instrument.
The instrument displays always the average of the last 6 results, in case the Damping is set to "6",
This makes readings more stable and protects for measuring outliers.
Damping can be set between "1" and "25".

Adjust Zero Offset

Adjust Zero Offset allows to set the instrument readings to "0".
Zero offset procedure requires clear almost turbidity free liquid inside of sensors flow cell.



Take Notice:

Use ultra pure liquids without bubbles / particles only in case you like to shift "0" reading of the instrument.

Adjust Gain

Adjust Gain allows a correction of the instrument reading at a known turbidity.

Example:

Instrument reading: "0.95"

Turbidity: "1.00"

Gain must be set to "1.05" to result a reading of "1.00" ($0.95 \times 1.05 = 1$)

Set Relay Setpoint

The "Relay Setpoint" defines the switching point of the alarm relay.

Set Relay Hysteresis

The hysteresis affects the relay setpoint.

Following settings: setpoint "2.00" and hysteresis "0.1"

Alarm will be switched at "2.1" and reset will be switched at "1.9"

So alarm switching point will be (setpoint + hysteresis)

and reset will be (setpoint – hysteresis).



Take Notice:

All Settings you change (Range, Damping, Zero, Gain, etc.) will be immediately transferred via USB to instrument's measuring electronic! But they are not in the permanent memory. So a reset or power on off will restore the old settings.

Apply Settings to Transmitter

Stores new settings (Range, Damping, Zero, Gain , etc.) to the permanent memory of the measuring electronic!

Save Settings to File

Store a backup of your instrument settings onto the hard disc of your computer / laptop. The backup files are in txt format and located at C:/Chemtronic_C1.



Take Notice:

We strictly recommend to backup the factory settings of the instrument before you do any changes at the instrument settings.

Read Settings from File

Read a backup file from the hard disc of your computer (data recovery).

Calibration



Danger, take notice!

Make sure that your pipe is not pressurized or at low pressure before you retract your probe from process. Before beginning any work the sensor must be cleaned and flushed thoroughly. Depending on customers' application product residues can be very dangerous (aggressive, poisonous).

Wear protective equipment if necessary (masks, safety glasses, protective gloves, etc.).

Notice:

- Qualified technical personnel must perform repair and maintenance.
- Avoid pulling and twisting of the probe cable.
- Avoid applying excessive force during installation and reinstallation of the probe.
- Please be sure that no dust or other particles penetrates into the optical assembly.
- In case of any dirty components please clean carefully by using fresh water and a cloth (lint-free), dry all components by using instrument air.
- Use suitable tools only.


Example calibration (Range 0-1000ppm)

- 1.) Fill a suitable vessel with low concentration calibration standard (zero liquid).
Typically distilled- or demineralized- water for the zeroing of turbidity probes.
Take notice that distilled- / demineralized water does not have an absolute zero concentration.
The turbidity of laboratory grade distillate shows typically a turbidity value of 0.1 - 0.3ppm
- 2.) Insert your properly cleaned probe into the Vessel.
Make sure that the optical path of the sensor is completely submerged.
Make sure that there are no bubbles in the optical path.
Protect inserted probe against ambient light.
- 3.) Connect your PC/Laptop via USB to the C1 Transmitter and start the C1-Config software.
- 4.) Use "Adjust Zero Offset" to set the instrument reading to a reading of "0ppm to 0.3ppm".
- 5.) Fill a suitable vessel with a 900ppm calibration standard.
Take notice, the calibration standard should have a concentration of about 75% to 100% of calibration range (0-1000ppm in our example).
- 6.) Insert your properly cleaned probe into the Vessel.
Make sure that the optical path of the sensor is completely submerged.
Make sure that there are no bubbles in the optical path.
Protect inserted probe against ambient light.
- 7.) Use "Adjust Gain" to set the instrument reading to a reading of "900ppm".

**Take Notice:**

Push button "Transfer Settings to Transmitter" in menu "Basic Settings" to finalize your calibration otherwise you will get the previous calibration after reset or power on/off. Transfer Settings to Transmitter will store the new calibration into the permanent memory of the instrument!

Measuring Screen

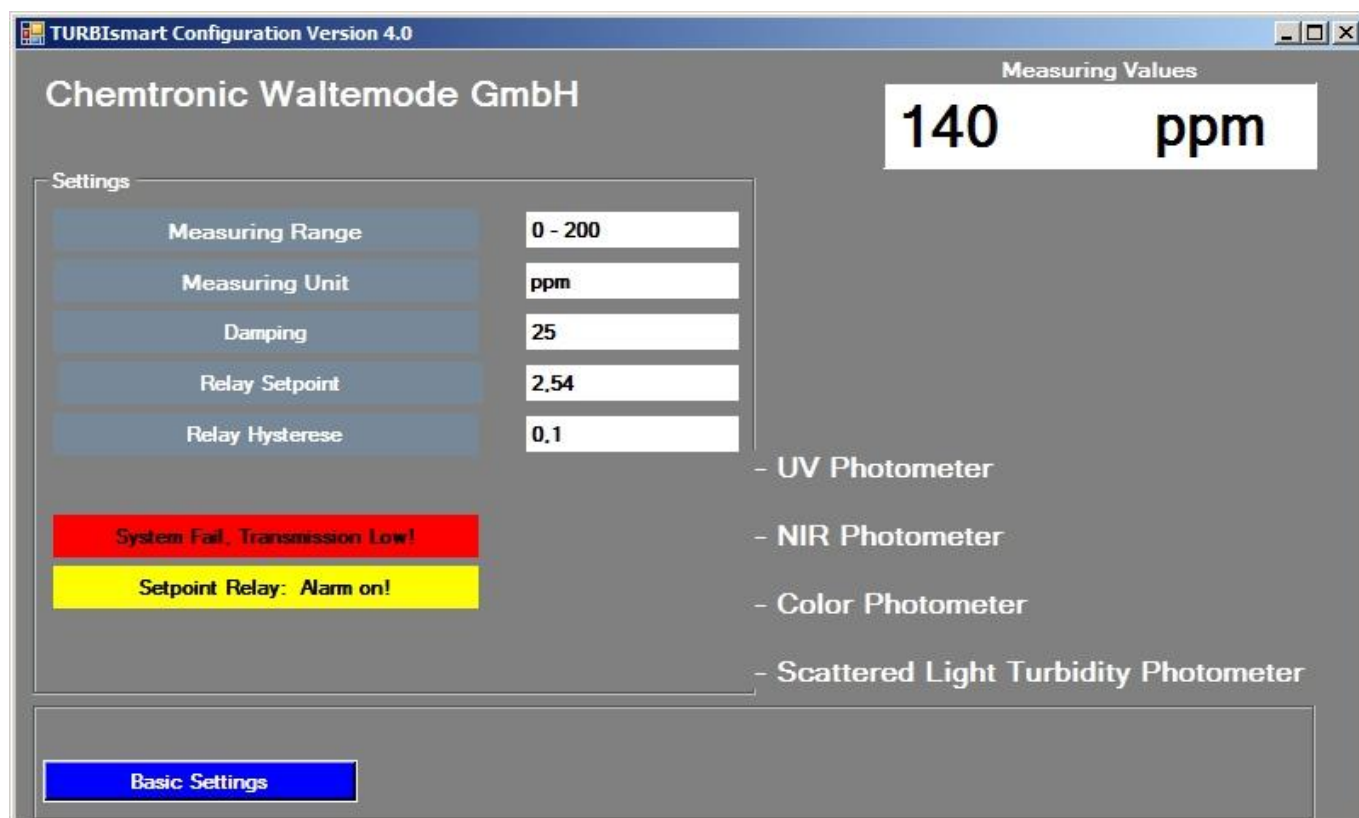

The measuring screen is the standard menu if the system is used with a PC / Panel PC. The menu displays the some general settings, system health information and measuring values.

Maintenance Work

The measuring transmitter model C1 is free of maintenance!

The probes models AP2, TURBImat-LC, TURBImat-HC have a low maintenance requirement which is limited on cleaning of the measuring windows.

Depending by application, coatings may form on the windows.

This may make it necessary to clean the measuring windows in fixed intervals.

The interval of the cleaning is determined by the operator depending on the application.

For cleaning, uninstall probe from process.

Clean the windows with a suitable cleaning agent, a cloth or a soft brush.

In case of mineral / lime deposits it is recommend clean the windows using diluted acetic acid.



Please take notice:

A factory inspection is recommended after 5 years of operation.



Please take notice:

In case of measuring hazardous substances (corrosive / toxic substances).

Appropriate protective measures must be taken before cleaning.

(Protective clothing, safety goggles, gloves, if necessary breathing mask).

Make absolutely sure that the measuring cell is emptied and depressurized before starting the cleaning procedure.

List of Spare Parts / Wearing Parts

No spare parts or wearing parts recommended.

Trouble Shooting

Complete failure: - Display and LEDs completely without function.

Remedy:

Check whether the supply voltage is present.

Use the wiring diagram to check whether the supply voltage is connected correctly.

For 24V DC operation, check polarity.

Check fuse F1 and replace it if necessary.

Contact Service.


Error: Check fuse - Fuse F2 for the measuring lamp defective.

Remedy:

Check that there is no short circuit in the lamp or the supply line to the lamp.

Check fuse F2 and replace it if necessary.

Contact service.

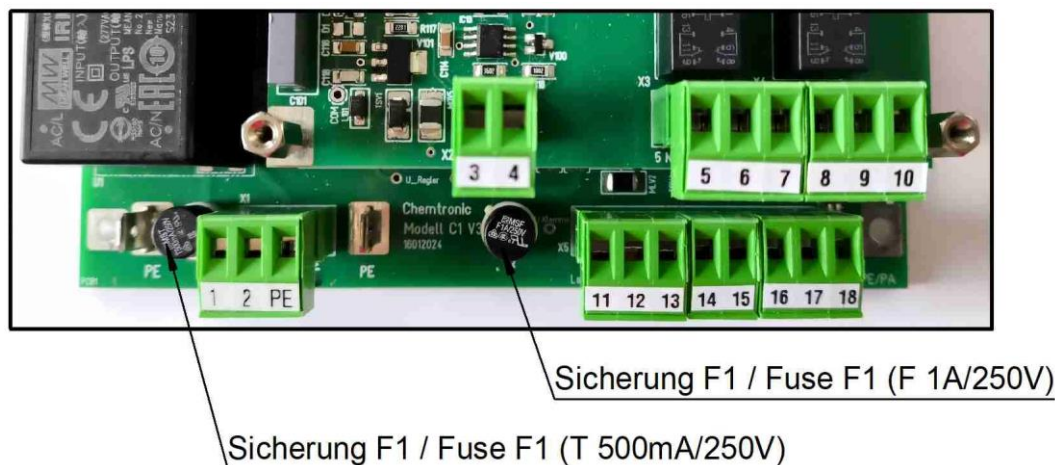

Error: Sensor LED - LED in the connected sensor is defective.

Error: Sensor lamp - Measuring lamp in the connected sensor is defective.

Remedy:

Use the wiring diagram to check whether the LED/measuring lamp is connected correctly.

Pay attention to polarity in case a LED- light source is connected.

Replace the LED/Lamp. (If necessary, follow the instructions in the relevant manual).

Contact service.

Error: Transmission Low - Sensor outside the physical measuring range, the light from the Led / measuring lamp in the sensor is completely absorbed.


Remedy: Abnormal process condition, solid concentration in the measured product extremely high.

Critical Absorption - Light absorption exceeds application-specific limit value.

Remedy:

Abnormal process condition, extremely high color/solids concentration in the sensor.

Terminal PC-Board Model C1

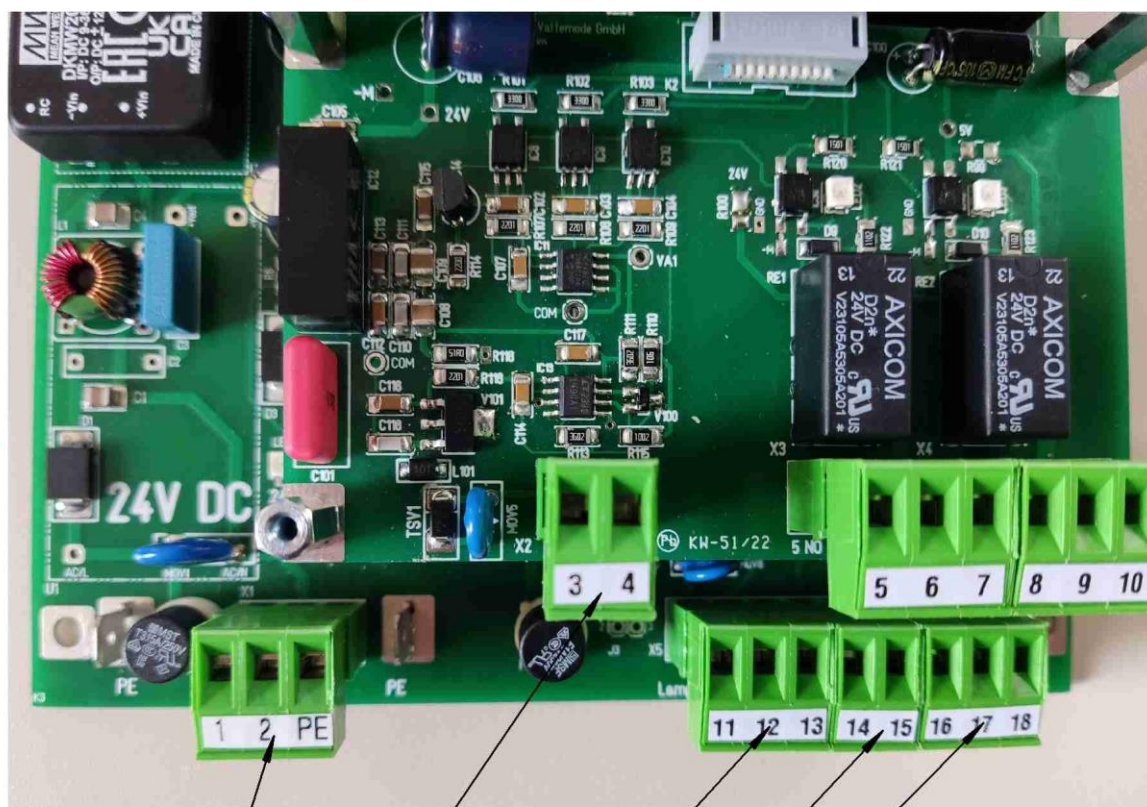

Supply Voltage
Versorgungsspannung


Analog out (0/4-20mA)
Analogausg. (0/4-20mA)

Lamp / LED
Lampe / LED

Signal 1 In


Signal 2 In

Dimensional Drawing Model C1

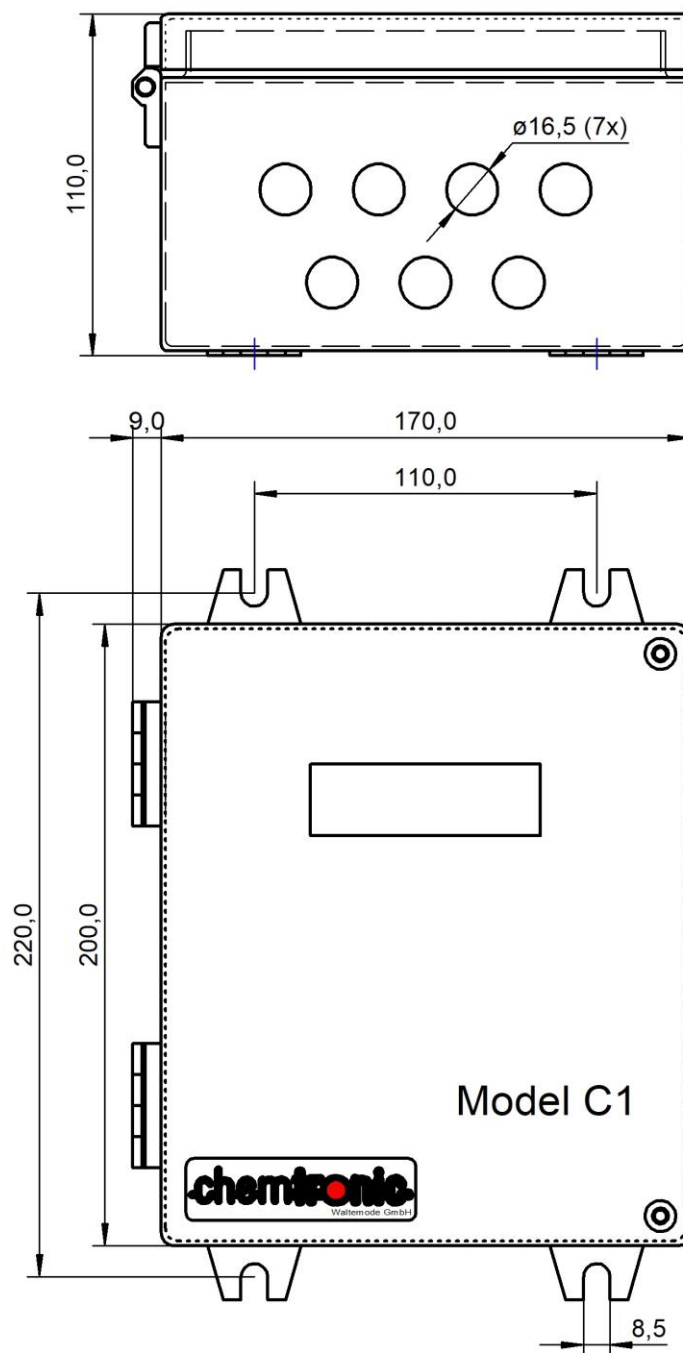


						Maßstab 1:2		Gewicht	
						Werkstoff 1.4404			
					Datum	Name		Abmessungen Modell C1	
				Bearb	07.12.2020	JHeidkamp			
				Gepr.					
				Norm					


Dimensional Drawings Model AP2

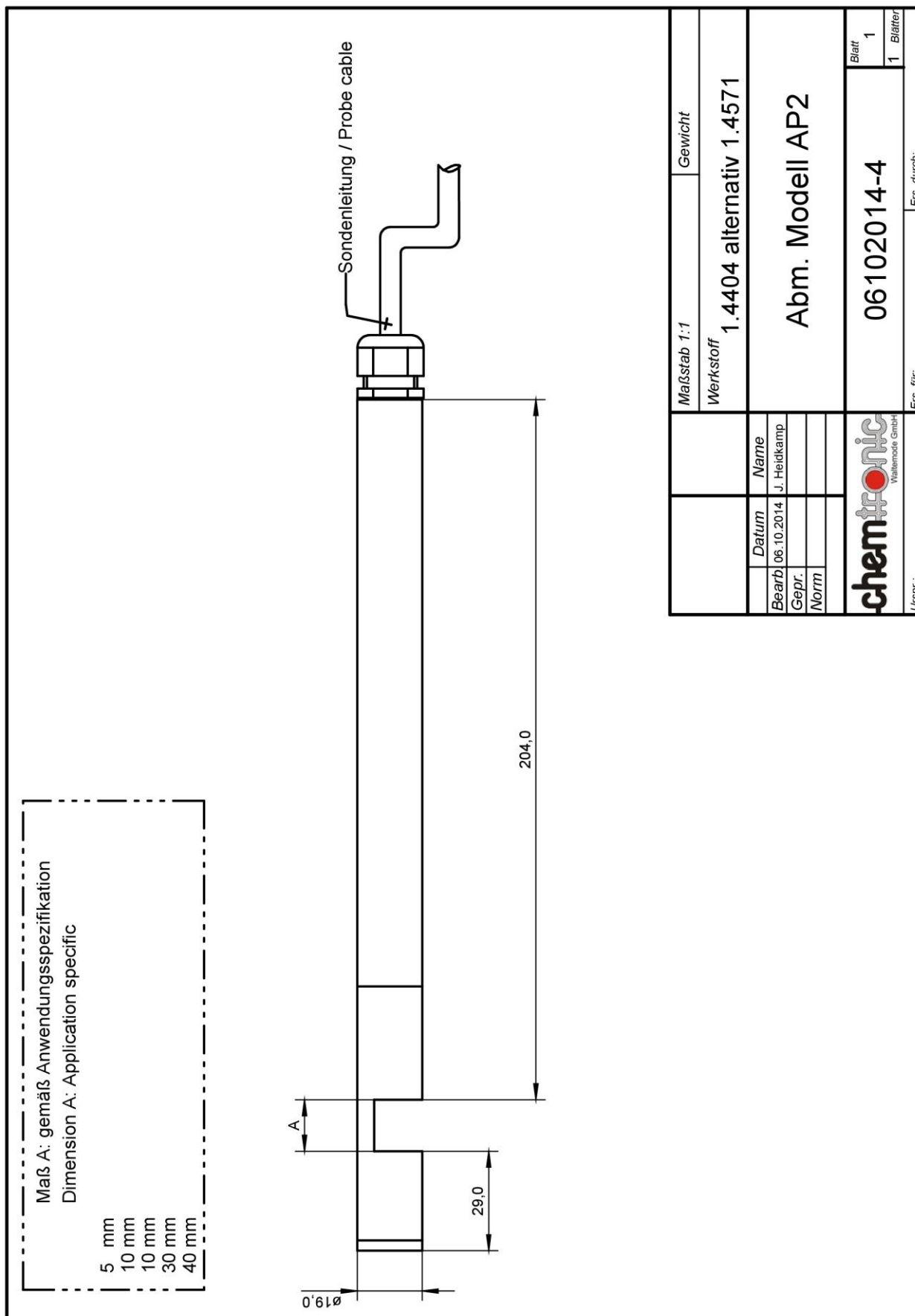


Dimensional Drawings Model AP2 with Flange

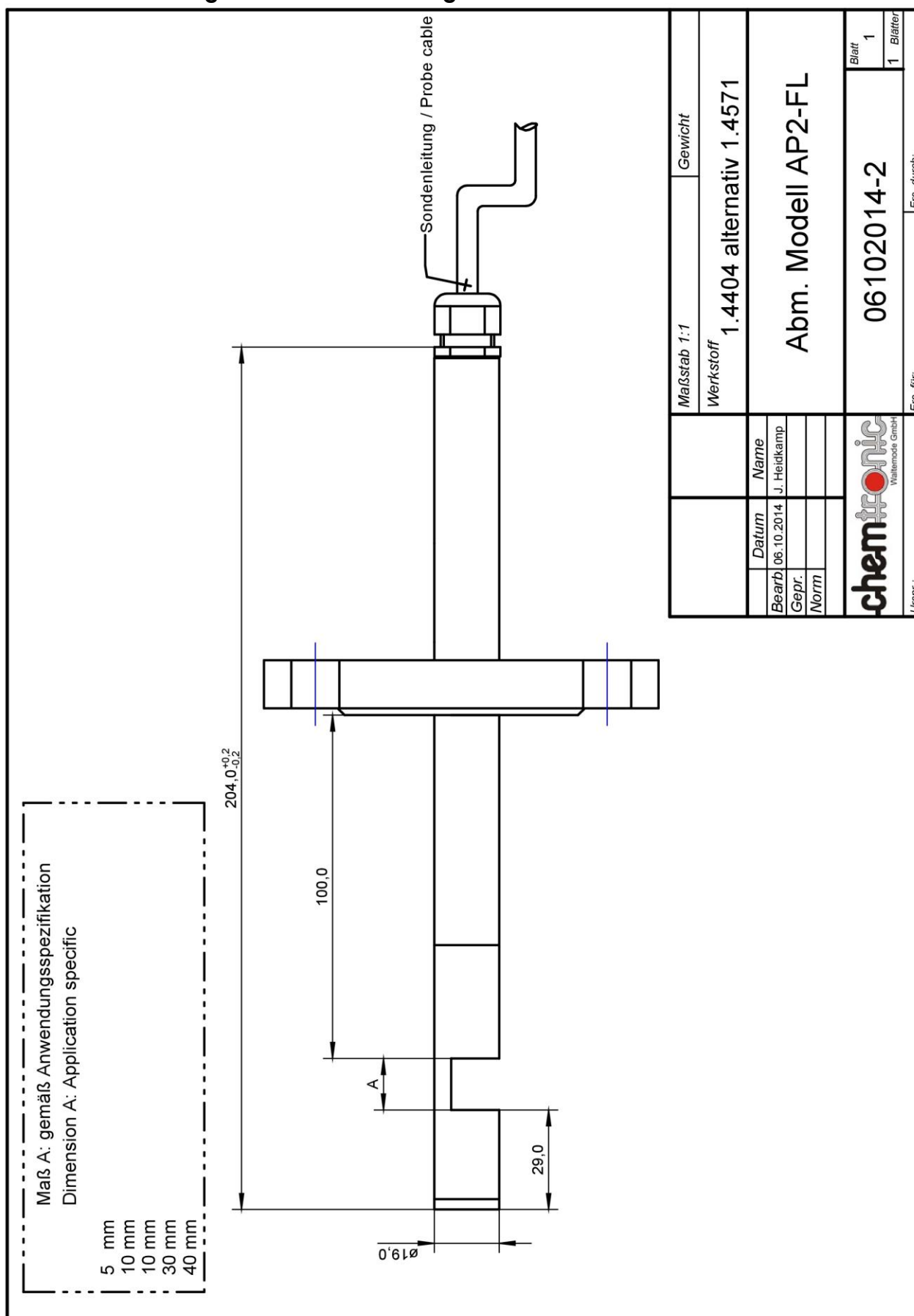

Dimensional Drawings Model TURBImat-HC

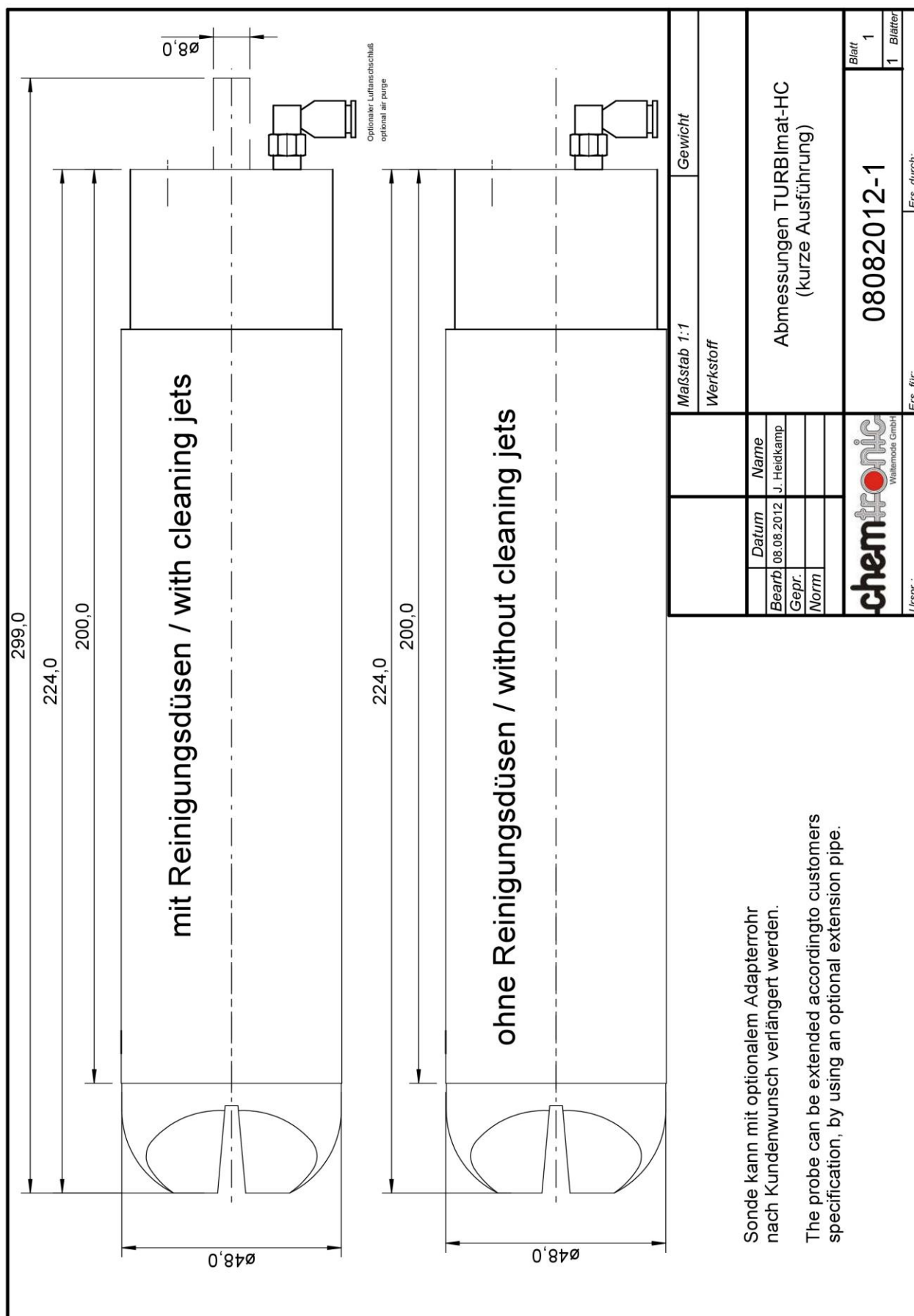


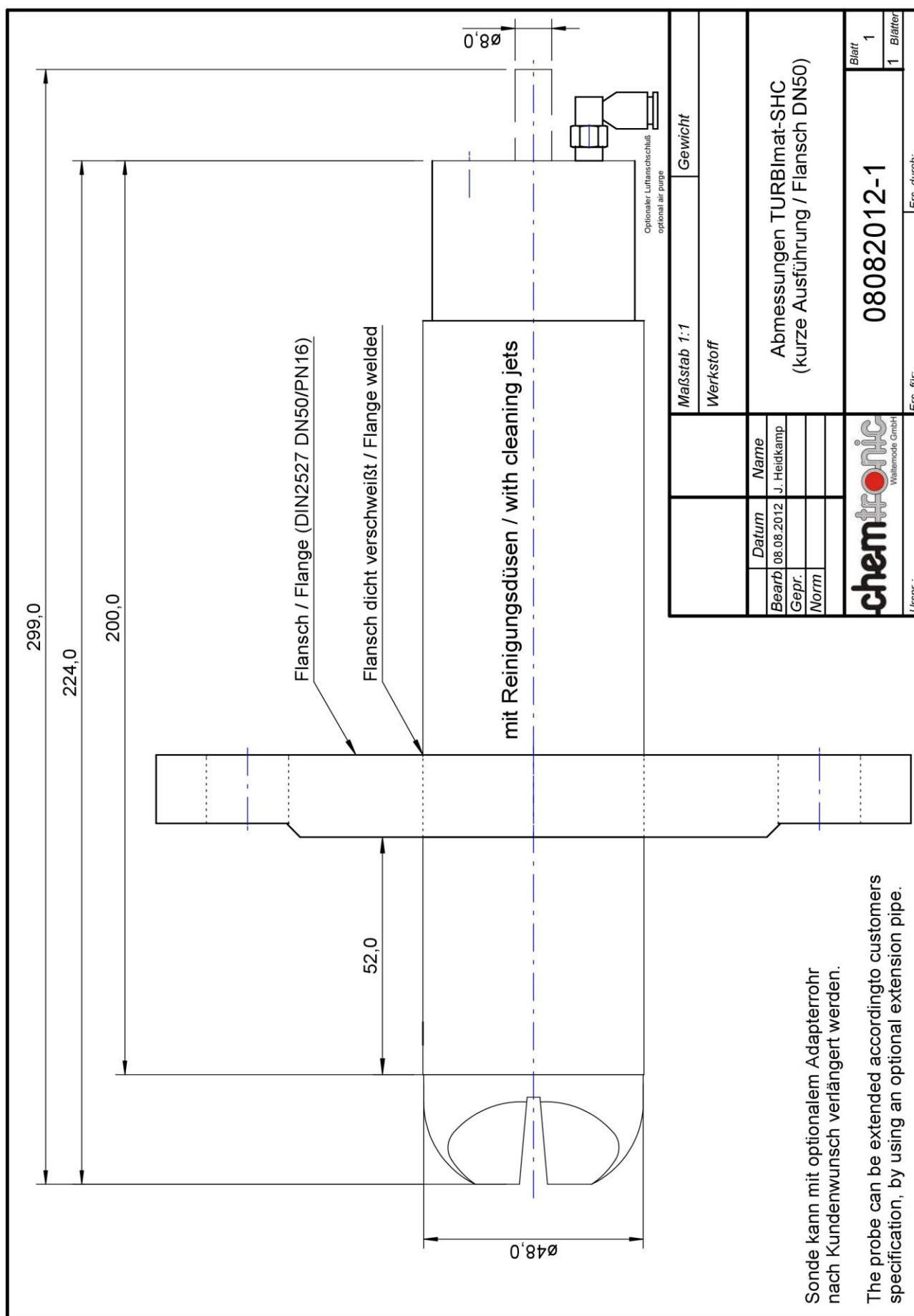

Dimensional Drawings Model TURBlmat-LC

Optionale Schaftlängen (L) / Optional shaft length (L)

150mm
700mm
1200mm
2000mm

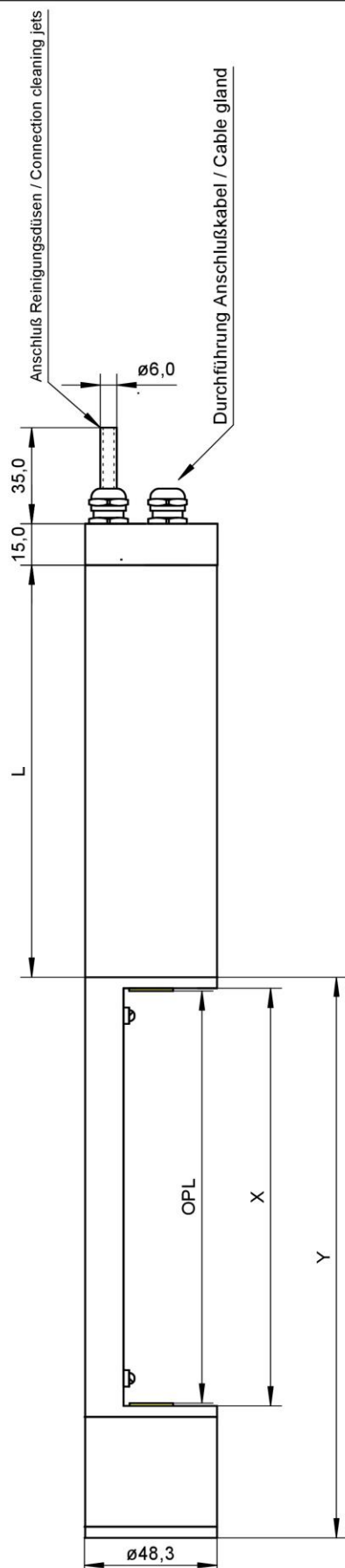



Optionale Schichtdicken (OPL) / Optional path length (OPL)

OPL	X	Y
10mm	12mm	64mm
30mm	32mm	84mm
50mm	52mm	104mm
100mm	102mm	154mm
150mm	152mm	204mm

Maßstab 1:2		Gewicht	
Werkstoff		1.4571 / 1.4404	
Abmessungen Modell TURBlmat-LC		Dimensional Drawing Model TURBlmat-LC	
chemtronic Waltemode GmbH		25012016-2	
Urspr.:		Ers. für:	
Datum		Name	
Bearb. 25.01.2016		J. Heidkamp	
Gepr.		Norm	
Blatt		1	
Blätter		1	

Data Sheet Transmitter Model C1

Transmitter Model C1

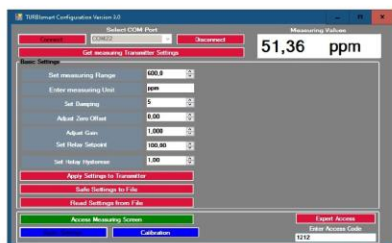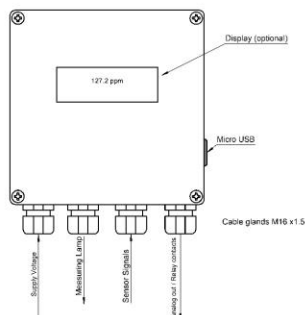

- Configuration via PC, Laptop
- Menu - based, intuitive User Interface
- Serial Interface: USB
- Optional backlit LCD
- Outputs: 4-20mA, Setpoint Relay, System fail Relay
- Multiple sample curve-fit
- Programmable Measuring Units: EBC, FTU, TEF, ppm
- Full backup of Configuration & Calibration

- Stainless Steel Enclosure (304SS / 1.4301)
- Optional Hazardous area Zone II (Ex nA IIA T3 Gc)*
- Optional Hazardous area Zone I alu enclosure (Ex d IIC T5/6)
- Protection class Enclosures: Nema 4x / IP65
- Stable and reliable excellent process control
- Suitable for all potable, waste and process waters

* Transmitter for hazardous area Zone II without Relays and external USB (non sparking device).

Description:

The transmitter model C1 can be used with all optical sensors of the Monitek / Chemtronic series. The programming/calibration of the system is done via a PC or notebook using the menu-driven software. Calibration is done with multiple samples. All settings and calibration values can be backed up, extern, at your PC's hard disc.

Applications:

- Scatter light turbidity measurement
- Absorption turbidity measurement
- Single channel color measurement
- Dual channel color measurement

Operational areas:

- Chemical industry
- Petrochemical industry
- Pulp & Paper
- Beer and beverages

Technical Data:

Temperature range: maximal 45°C

Enclosure material: 1.4301

Measuring range: Sensor specific (programmable)

Interface: USB

Software: Configuration/Calibration via Windows PC

Reproducibility: $\pm 1 \%$

Protection class: IP65 / NEMA 4X

Supply: 18-27 VDC optional 90-260VAC / 50/60Hz

Analog output: 4-20mA, 500 Ohm (isolated)

Relay outputs: Contact ratings 24VDC/1A

Data Sheet Probe Model AP2

Turbidity Probe Model AP2



- 19 mm probe technology
- Measurement method: NIR- Absorbance
- Different measurement wavelength on request
- Low maintenance
- Extended calibration interval: Typical 12 month
- Sapphire measuring windows
- Hygienic design
- Cleaning in place (CIP)
- Process installation via TH- Varivent plates available
- Process installation via DIN / Ansi flange available
- Process insertion via ball valve with adapter available

Description:

The probe model AP2 uses the principle of near infrared light absorption to detect turbidity in liquids. The basic version of the probe allows for process temperatures of max. 45°C. The optional high temperature version allows for 105°C permanent temperatures. The transmitter model C, C1 or CCAb is required to use the AP2 probe. The system has been designed for continuous operation with long life time. The 19mm probe design allows insertion to laboratory fermentation tanks. Different flanges, TH- Varivent adaption, Tri Clover adaption or ball valve installation are available to allow application specific insertion of the probe in nearly any process pipe.

Applications:

- Product Separation
- Turbidity measurement
- TSS concentration measurement
-

Operational areas:


- Breweries
- Beverages
- Dairies
- Bio technology
- Water / Waste water
- Chemical industry

Technical Data:

Line size:	DN 25 / 1" up to	Measurement ranges:	typical 0–100EBC, 0–500EBC
Process pressure:	PN16 / ANSI class 300	Reproducibility:	± 1 %
Process temperature:	45°C (131°C short time CIP)	Measurement wavelength:	880nm (other on request)
Optional:	105°C permanent temperature	Protection:	IP65 / NEMA4X
Sensor material:	1.4404 / 316L	Hazardous Area:	ATEX, Zone 1 & Zone 2
Sight glass material:	Sapphire	Sterilization:	CIP (optional SIP)
Gasket material:	EPDM (other on request)		

Data Sheet Probe Model TURBImat-HC

Model TURBImat-HC (TBMHC)

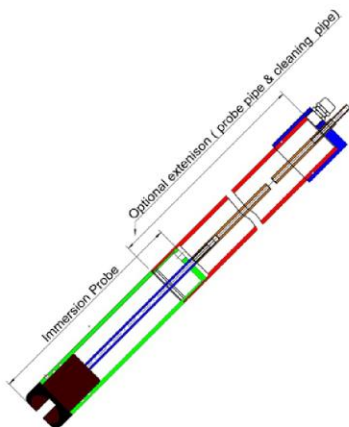

- Immersion probe
- Low maintenance
- Extended calibration interval: Typical 6 month
- Multiple probe extensions available
- Immersion / or pipe installation via 2" insertion adapter
- Insertion / removal without process interruption
- Multiple optical path lengths available (max. 10 mm)
- Optional available with DIN- or ANSI flanges
- Optional 316 SS measurement head / sapphire windows
- Optional cleaning jets
- Wide span of range

Description:

The turbidity probe model TURBImat-HC uses the principle of light absorption to detect dissolved or suspended particles in liquids. The basic probe is 230mm long and will be immersed into the measured product.

Extension tubes allow the installation to pipes as well as the immersion in tanks or open channels. Pipe installation is typically done by using a 2" ball valve or a mounting flange at the probe shaft.

The sensor requires the transmitter model C, C1 or CCAb and is suitable for the measuring of high concentrations in liquids. The calibration depends by application, measuring range and unit.

Applications:

- Product concentration
- Sedimentation
- Flotation
- Sludge concentration

Operational areas:


- Chemical industry
- Petrochemical industry
- Pulp & Paper
- Power plants & Waste water

Technical Data:

Line size:	minimum diameter DN 80	Immersion depth:	maximum 80m
Process pressure:	maximum 8Bar	Measurement range:	typical 0-6%
Process temperature:	maximum 60°C	Reproducibility:	± 1 %
Probe material:	1.4471 / 316SS	Detector:	Silica diode
Measurement head:	PVC (optional 1.4571 / Sapphire)	Measurement wavelength:	typical 880nm
Gasket material:	Application specific	Protection:	IP68
Cleaning:	Optional cleaning jets (air/water)	optional hazardous area:	ATEX Zone I / Zone II

Data Sheet Probe Model TURBImat-LC

Model TURBImat-LC (TBMLC)

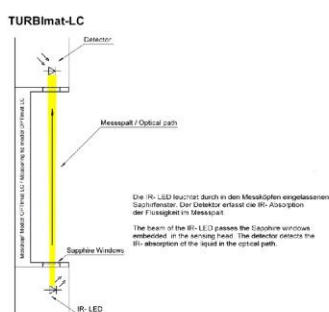

- Low maintenance
- 30cm immersion probe
- Pipe extension 0,5m for pipe insertion via ball valve
- Insertion/removal without process interruption
- Pipe extension up to 2m for installation in open channel
- Calibration interval: Typical 6 month
- Wide span of measurement ranges
- Sight glass material: Sapphire
- Sight glass cleaning via cleaning jets
- Optional available with DIN- or ANSI flanges

Description:

The turbidity probe model TURBImat-LC uses the principle of light absorption to detect dissolved or suspended particles in liquids. The transmitter model C, C1 or CCAB is required to use this probe. The system has been designed for continuous operation with long life time. The probe design allows installation into large process pipes, tanks or open channels. Up to a process pressure of 3bar and non critical product, process insertion and removal of the probe can be easily done via a 2"ball valve. So calibration and maintenance ensues without process interruption. Other process connections like flanges, swivel clamps, etc. are available on request.

Applications:

- Product concentration
- Sedimentation
- Flotation
- Oil in water / Water in oil

Operational areas:

- Chemical industry
- Petrochemical industry
- Pulp & Paper
- Power plants

Technical Data:

Line size:	minimum diameter DN 200	Measurement range:	0–100ppm, 0–2000ppm
Process pressure:	maximum PN 16 / ANSI class 150	Reproducibility:	± 1 %
Process temperature:	maximum 60°C	Detector:	Silica diode
Probe material:	1.4471 / 316SS	Measurement wavelength:	typical 880nm
Sight glass material:	Sapphire	Protection:	IP65 / NEMA 4X
Gasket material:	Application specific	optional hazardous area:	ATEX Zone I / Zone II

EU Declaration of Conformity (CE)

EU-Konformitätserklärung EU Declaration of Conformity

Der Hersteller:
The Manufacturer:
Chemtronic Waltemode GmbH
Niederstraße 14
40789 Monheim am Rhein



erklärt hiermit das das Prozess Fotometer
herewith declares that the process photometer:

Messverstärker / Transmitter Mod.:	C1 in Verbindung mit / with
Sensoren Modelle / Sensors & Probes Mod.:	AP2, CSW, LAS, MTC., MSC., TBM,

folgenden europäischen Richtlinien entspricht:
complies with the following European directives:

2014/30/EU	Richtlinie über die elektromagnetische Verträglichkeit Electromagnetic compatibility
2014/35/EU	Niederspannungsrichtlinie Electrical equipment designed for use within certain voltage limits

Angewandte harmonisierte Normen:
Applied harmonized standards:

EN61010-1 /2020	Sicherheitsbestimmungen für Mess-, Steuer-, Regel- und Laborgeräte Safety requirements for equipment for measurement, control, & laboratory use
EN61000-6,2 /2019	Störfestigkeit für Industriebereiche Immunity standard for industrial environments
EN61000-6,4 /2011	Elektromagnetische Verträglichkeit (Störaussendung) Emission standard for industrial environments

Monheim den 16.06. 2020

Gunnar Waltemode
Geschäftsführer / General Manager